

# Magnetically actuated glaucoma drainage device with adjustable flow properties after implantation

**Citation for published version (APA):**

Pereira, I. C. F., Wyss, H. M., Beckers, H. J. M., & den Toonder, J. M. J. (2021). Magnetically actuated glaucoma drainage device with adjustable flow properties after implantation. *Engineering Proceedings*, 4(3). <https://doi.org/10.3390/Micromachines2021-09584>

**DOI:**

[10.3390/Micromachines2021-09584](https://doi.org/10.3390/Micromachines2021-09584)

**Document status and date:**

Published: 16/04/2021

**Document Version:**

Publisher's PDF, also known as Version of Record (includes final page, issue and volume numbers)

**Please check the document version of this publication:**

- A submitted manuscript is the version of the article upon submission and before peer-review. There can be important differences between the submitted version and the official published version of record. People interested in the research are advised to contact the author for the final version of the publication, or visit the DOI to the publisher's website.
- The final author version and the galley proof are versions of the publication after peer review.
- The final published version features the final layout of the paper including the volume, issue and page numbers.

[Link to publication](#)

**General rights**

Copyright and moral rights for the publications made accessible in the public portal are retained by the authors and/or other copyright owners and it is a condition of accessing publications that users recognise and abide by the legal requirements associated with these rights.

- Users may download and print one copy of any publication from the public portal for the purpose of private study or research.
- You may not further distribute the material or use it for any profit-making activity or commercial gain
- You may freely distribute the URL identifying the publication in the public portal.

If the publication is distributed under the terms of Article 25fa of the Dutch Copyright Act, indicated by the "Taverne" license above, please follow below link for the End User Agreement:

[www.tue.nl/taverne](http://www.tue.nl/taverne)

**Take down policy**

If you believe that this document breaches copyright please contact us at:

[openaccess@tue.nl](mailto:openaccess@tue.nl)

providing details and we will investigate your claim.

ICMA  
2021

# 1st International Conference on Micromachines and Applications

15–30 APRIL 2021 | ONLINE



micromachines



## Magnetically actuated glaucoma drainage device with adjustable flow properties after implantation

Inês Pereira<sup>1,2,\*</sup>, Hans Wyss<sup>1,2</sup>, Henny Beckers<sup>3</sup>, and Jaap den Toonder<sup>1,2</sup>

<sup>1</sup> Microsystems Research Section, Eindhoven University of Technology, Eindhoven, The Netherlands

<sup>2</sup> Institute for Complex Molecular Systems, Eindhoven University of Technology, Eindhoven, The Netherlands

<sup>3</sup> Maastricht University Medical Center, University Eye Clinic, Maastricht, The Netherlands

\* Corresponding author: [i.c.figueiredo.pereira@tue.nl](mailto:i.c.figueiredo.pereira@tue.nl)


# Abstract

Glaucoma is a leading cause of preventable blindness worldwide. A rise in the intraocular pressure (IOP) is a major risk factor for this disease, and results from an elevated resistance to aqueous humor outflow from the anterior chamber of the eye. Glaucoma drainage devices provide an alternative pathway through which aqueous humor can effectively exit the eye, thereby lowering the IOP. However, post-operative IOP is unpredictable and current implants are deficient in maintaining IOP at optimal levels. To address this deficiency, we are developing an innovative, non-invasive magnetically actuated glaucoma implant with a hydrodynamic resistance that can be adjusted following surgery. This adjustment is achieved by integrating a magnetically actuated microvalve into the implant, which can open or close fluidic channels using an external magnetic stimulus. This microvalve was fabricated from poly(styrene-*block*-isobutylene-*block*-styrene), or 'SIBS', containing homogeneously dispersed magnetic microparticles. "Micro-pencil" valves of this material were fabricated using a combination of femtosecond laser machining with hot embossing. The glaucoma implant is comprised of a drainage tube and a housing element fabricated from two thermally-bonded SIBS layers with the microvalve positioned in between. Microfluidic experiments involving actuating the magnetic micro-pencil with a moving external magnet confirmed the valving function. A pressure difference of around 6 mmHg was achieved which is sufficient to overcome hypotony (i.e. too low IOP) – one of the most common post-operative complications following glaucoma surgery.


**Keywords:** Glaucoma; Intraocular pressure (IOP); Adjustable glaucoma implant; Magnetic microvalve; Microfabrication.

# Introduction – What is glaucoma?

- Eye disease characterized by optic nerve damage and visual field loss
- Leading cause of irreversible blindness worldwide, with over **70 million people affected!**
- Most important known risk factor: **elevated intraocular pressure (IOP)**

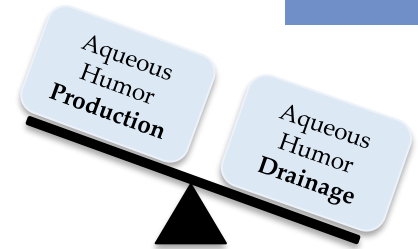


## Aqueous humor dynamics:

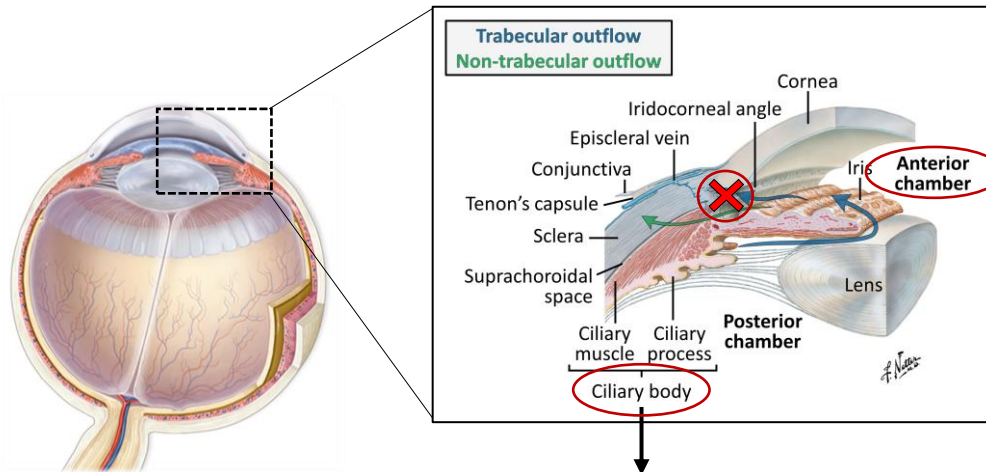

Continuous production  
of aqueous humor  
(2,5  $\mu$ L/min)

High IOP damages the  
optic nerve

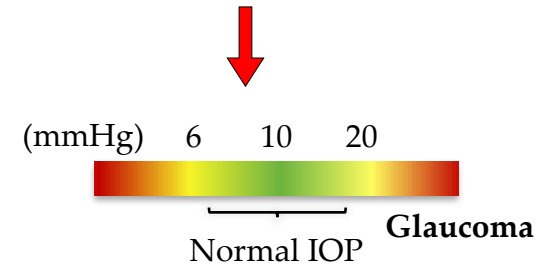


# Introduction – Current treatment options

- ✓ Pharmacological medication
  - ✓ Laser treatment
  - ✓ Surgery – ↑ Aqueous humor drainage
- } ↓ Aqueous humor production OR ↑ Aqueous humor drainage

**Glaucoma drainage devices**  
 Bypass the eye's natural outflow pathways to provide an alternative/artificial route for aqueous humor to effectively exit the eye, thereby reducing IOP


**Baerveldt® implant**

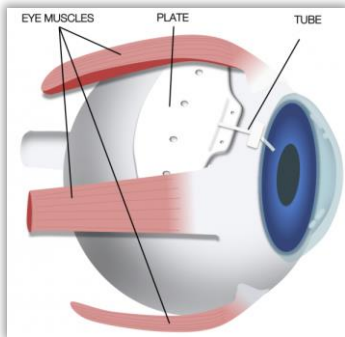

**The Ahmed® Glaucoma Valve**



**PRESERFLO™ MicroShunt**

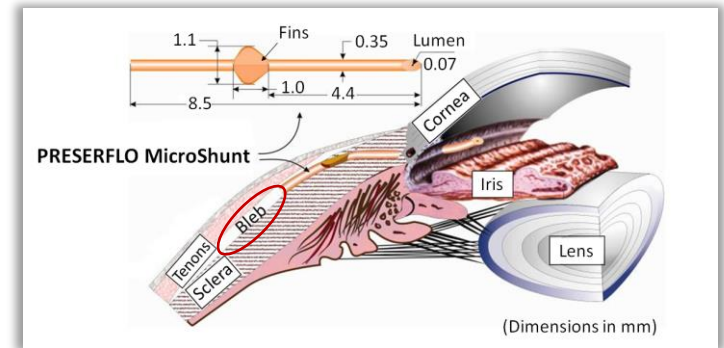


- ✓ Most glaucoma implants drain aqueous humor into the so-called “bleb”
- ✓ Aqueous humor in the bleb is then absorbed by the subconjunctival microvasculature

# Introduction – Glaucoma drainage devices

## Most common postoperative complications

- ✓ **Hypotony** – 5 mmHg or less (usually occurs in the early postoperative period)
- ✓ **Bleb encapsulation** (occurs weeks/months after surgery)

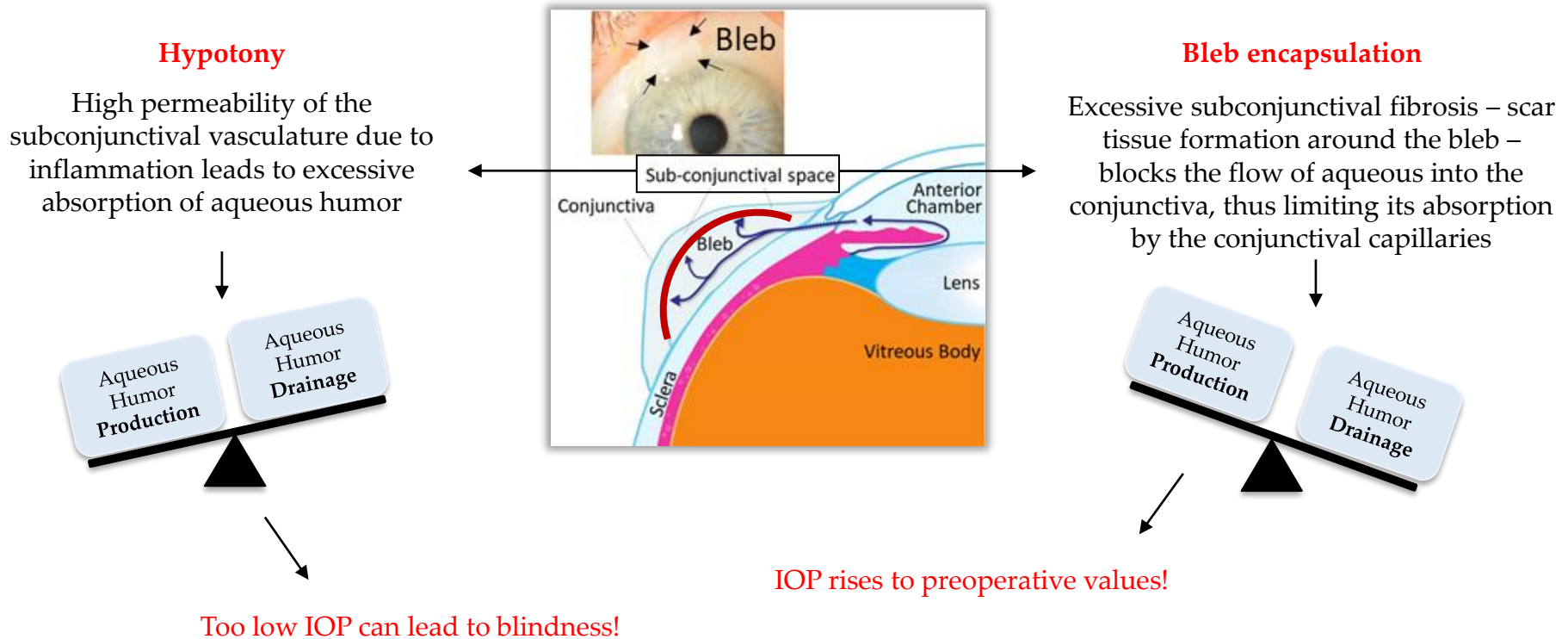



# Introduction – Aim of the project

- ✓ Current glaucoma implants provide always the same fluidic resistance; thus, fluid flow cannot be adapted according to each patient's needs (**not all patients suffer from hypotony or bleb encapsulation!**)

Solution



- ✓ Design and fabricate a novel **SMART GLAUCOMA DRAINAGE DEVICE** with **flow control** properties

Flow of aqueous humor through the implant will be adjusted if necessary, after surgery and non-invasively, through the actuation of a **stimuli-responsive valve** integrated in the implant!


- ✓ Fabricate the stimuli-responsive microvalve and develop smart solutions to integrate it into the implant
- ✓ Evaluate the flow control properties of the smart implant by carrying out microfluidic experiments
- ✓ **Final goal:** produce potential implant prototypes that can be used for animal testing

# Magnetically adjustable glaucoma implant

- ✓ Smart implant is comprised of a drainage tube and a housing element, both made from a highly biocompatible, bioinert material called **poly(styrene-*block*-isobutylene-*block*-styrene), or SIBS**:

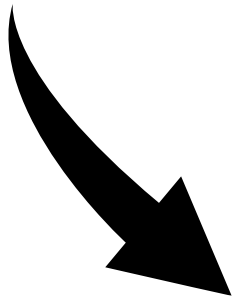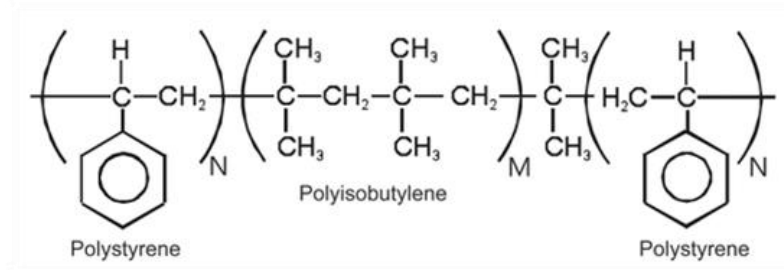




Smart implant:

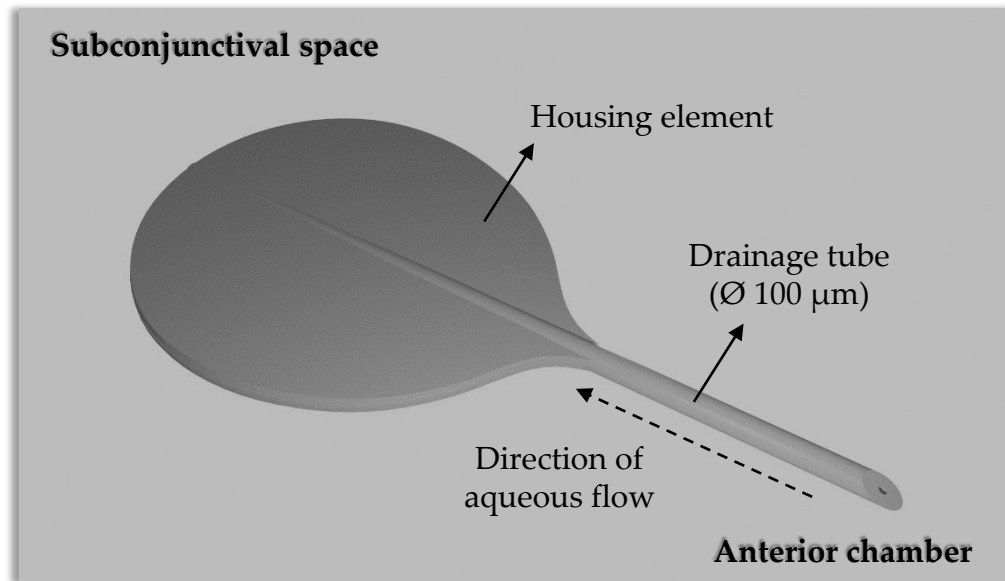




# Magnetically adjustable glaucoma implant

- ✓ The housing element of the smart implant contains an actuation chamber where the microvalve, with the shape of a “micro-pencil”, is positioned
- ✓ An external magnet is used to move this microvalve to the open or closed positions: (i) when open, the device is in a “high flow” mode; (ii) when closed, the device is in a “low flow” mode:

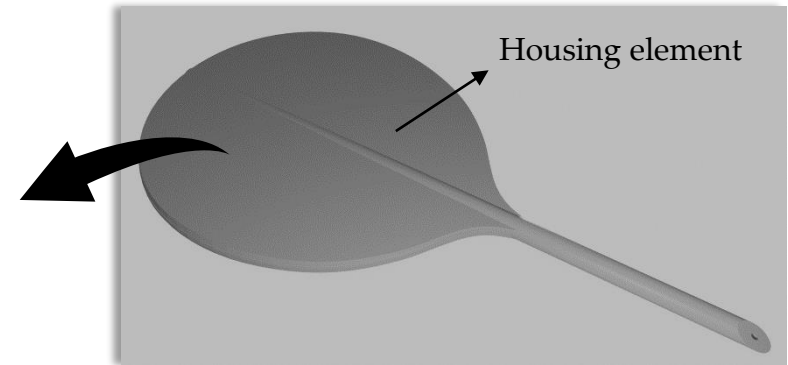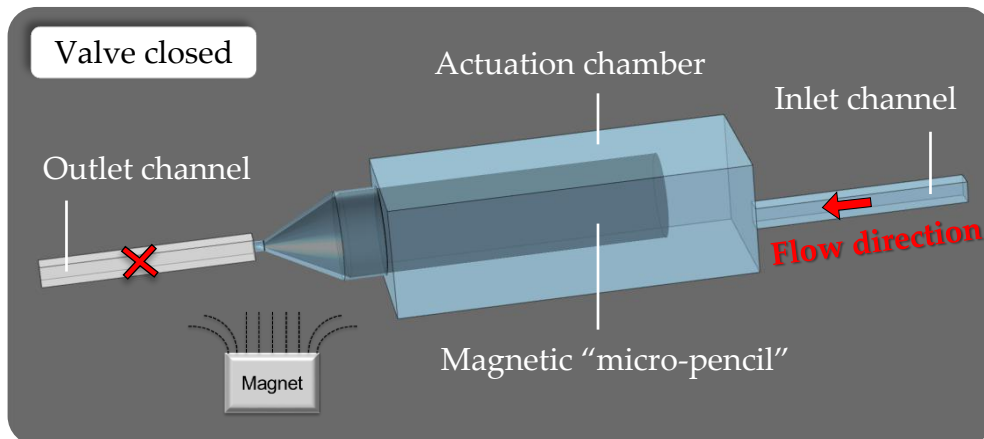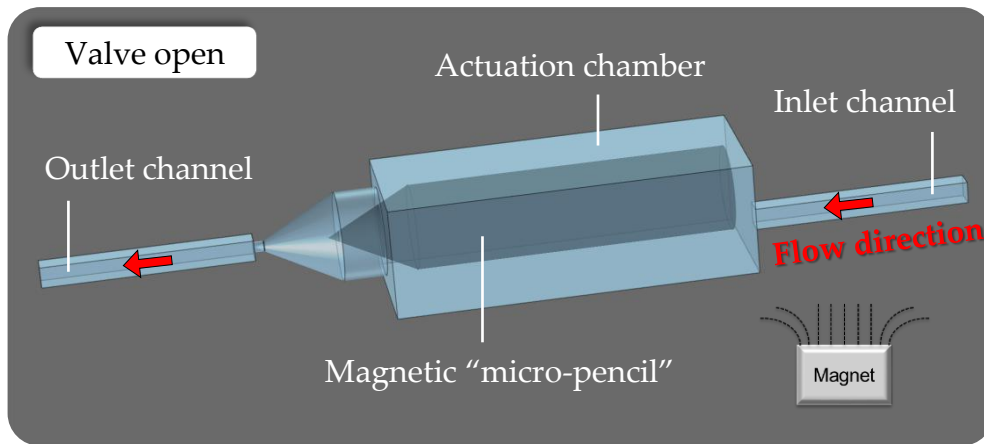

# Fabrication of the magnetic microvalve

- ✓ The magnetic microvalve is made from SIBS mixed with carbonyl iron powder (CIP): SIBS:CIP = 1:1 (w/w)
- ✓ Hot melt extrusion with mixing was used to fabricate the batch material (magnetic SIBS)

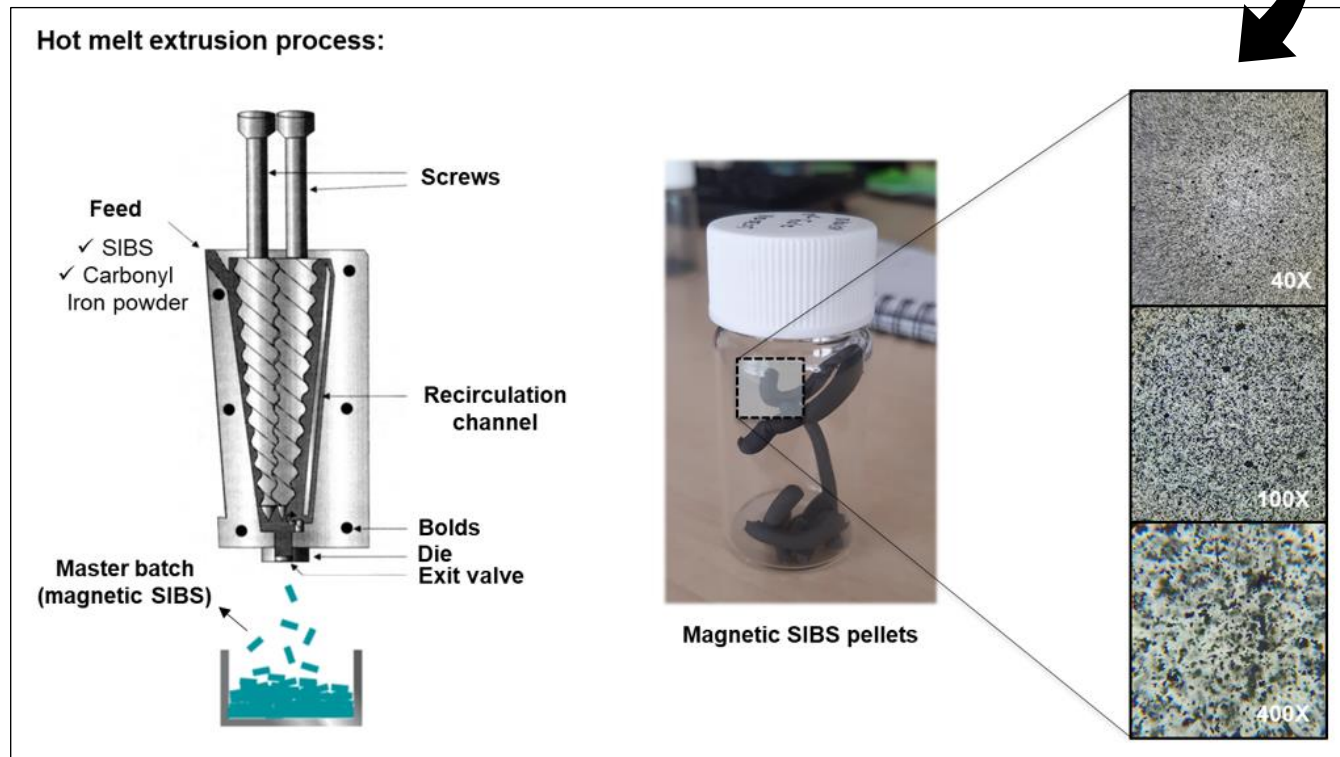


# Fabrication of the magnetic microvalve

- ✓ The micro-pencil valve mold was fabricated using the FEMTOprint




## Femtosecond laser machining process:

- ✓ Modifying a material by exposing to femtosecond laser in a desired pattern followed by its etching to complete the machining process.

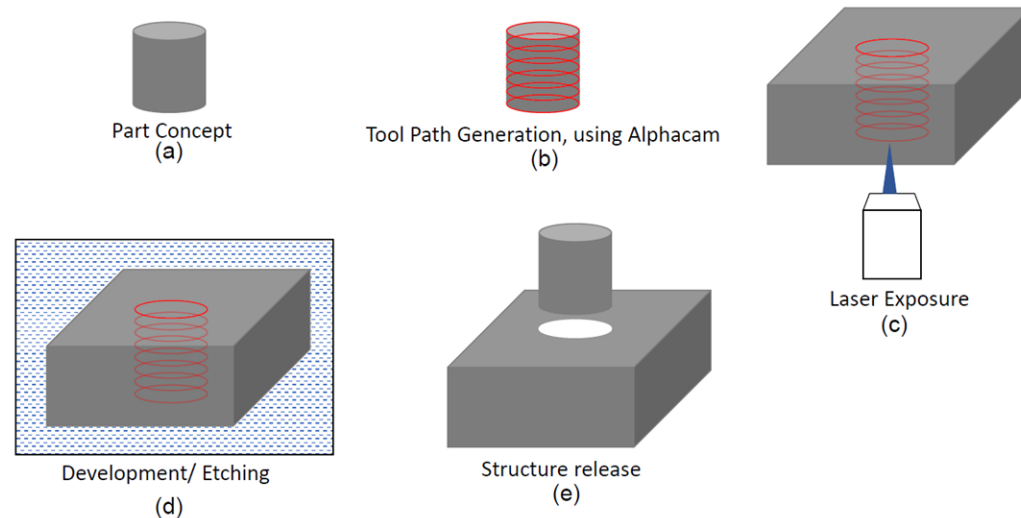

# Fabrication of the magnetic microvalve

- ✓ FEMTOprint micro-pencil mold

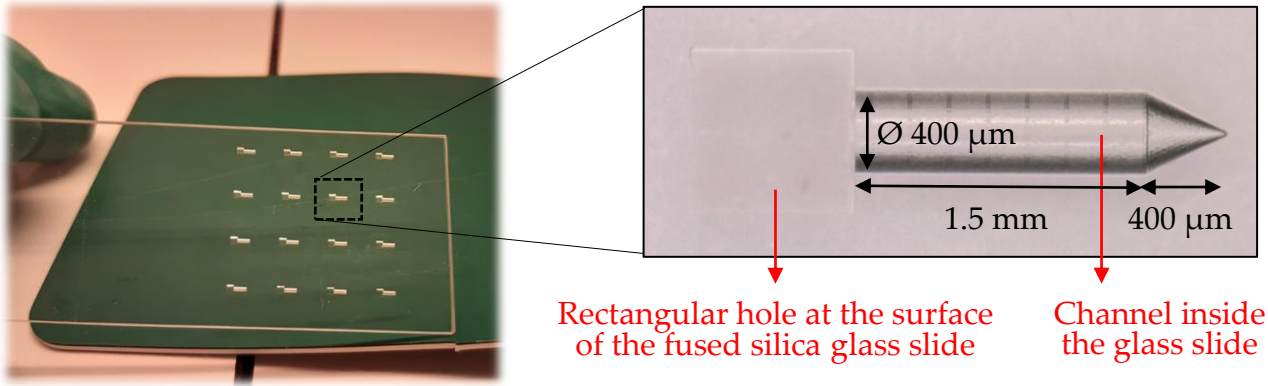

## Machining and etching times:

- ✓ Machining  $\approx$  6h
- ✓ KOH etching  $\approx$  48h

## Magnetic Micro-Pencil:




## Hot embossing:


- ✓ FEMTOprint micro-pencil mold and Magnetic SIBS pellets are “sandwiched” between two brass plates
- ✓ 150 °C, 15 min, 5 tons

# Fabrication of the housing element

- ✓ Housing element fabricated from two thermally-bonded SIBS layers with channels
- ✓ Molds for top and bottom layers were fabricated using the FEMTOprint

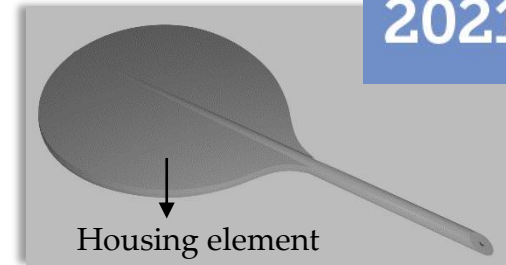

## FEMTOprint mold for top layer:

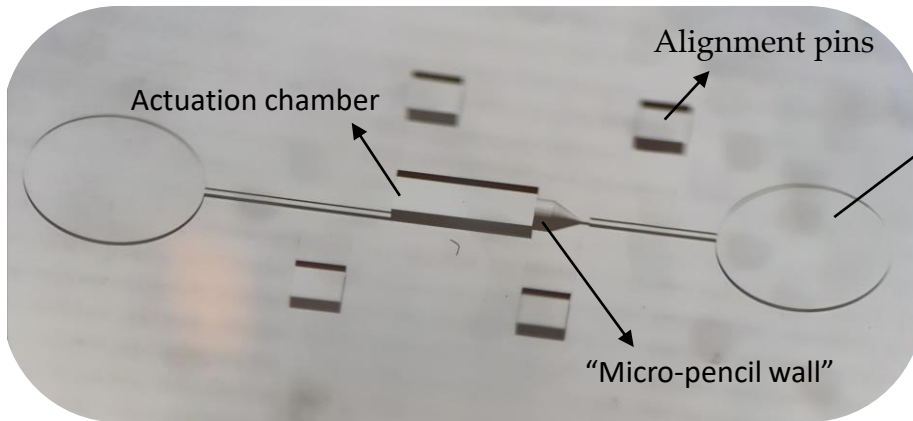


Inlet and outlet ports used for *in vitro* microfluidic experiments

## Machining and etching times:

- ✓ Machining  $\approx$  12h
- ✓ KOH etching  $\approx$  7h

## Replica-molded SIBS top layer:


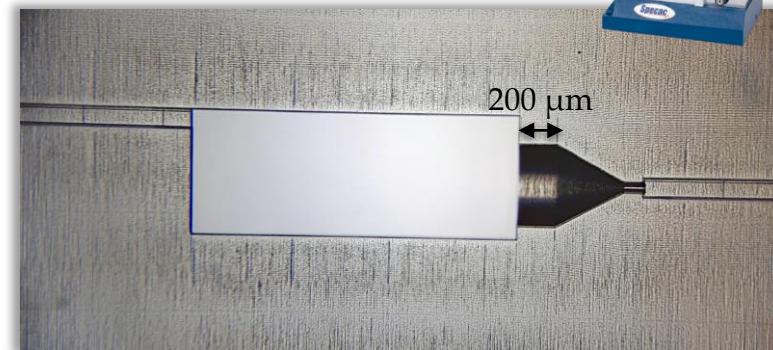


## Hot embossing:

- ✓ 150 °C, 15 min, 5 tons

# Setup for microfluidic experiments

- ✓ Aim: evaluate how the pressure upstream the implant (IOP) changes when the microvalve is moved to the closed or open positions.

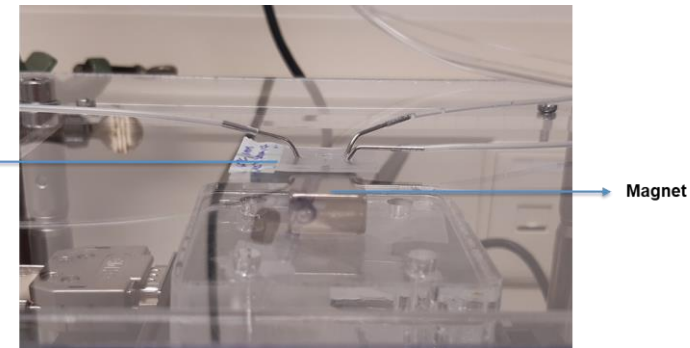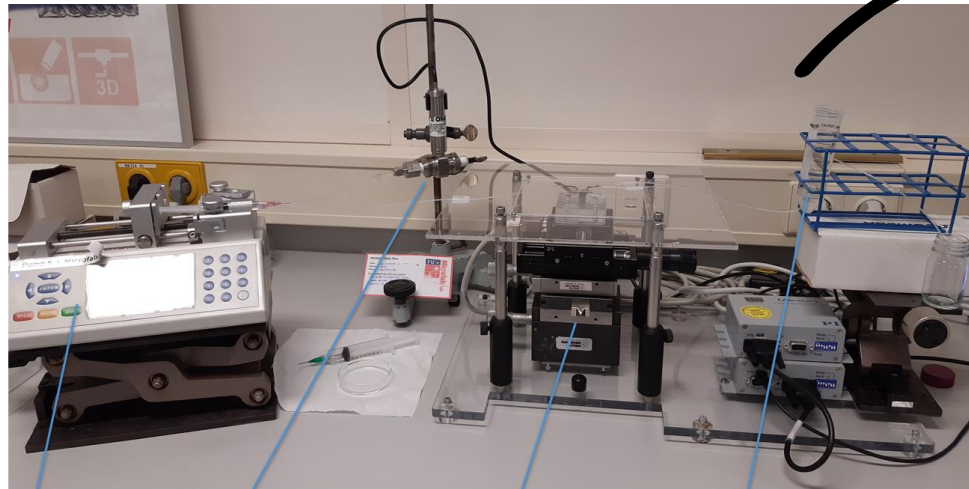

**Syringe pump:**  
pumping DI water  
at 2.5  $\mu\text{L}/\text{min}$

**Pressure transducer:**  
measuring the pressure  
upstream the device (inlet)


**X,Y and Z stages:**  
move the magnet in one direction  
at a pre-defined velocity

Outlet of the device connected to a column of  
water of 6.8 cm, which corresponds to 5  
mmHg (mimicking the pressure in the bleb)

It is expected that the pressure  
upstream the device increases when  
the microvalve is moved to the closed  
position!



# Microfluidic experiment result

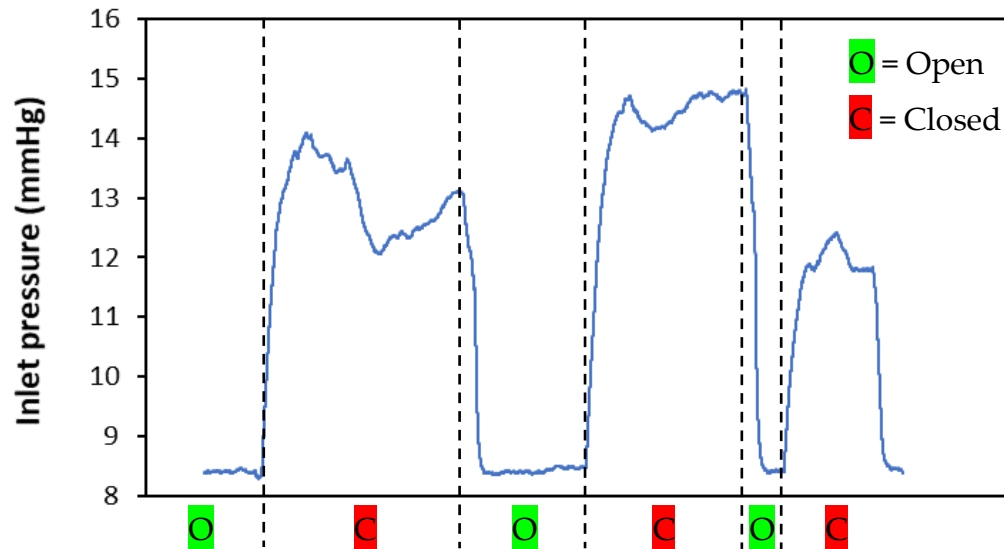

- Closed position: **IOP increase up to 6 mmHg**
- Pressure difference achieved is sufficient to overcome hypotony (i.e. too low IOP)
- Pressure difference varies slightly between open/closed cycles

# Conclusions

- ✓ The magnetic SIBS batch material was successfully fabricated using the hot melt extrusion process
- ✓ The femtosecond laser machining process is an effective technique to fabricate the molds for both the device and the microvalve
- ✓ Microfluidic tests have shown that, when in the closed state, the magnetically actuated microvalve can provide a sufficient hydrodynamic resistance that leads to an increase in the pressure upstream the implant (IOP)
- ✓ Pressure difference achieved is sufficient to overcome hypotony

# Future work

- ✓ Leaching tests, as well as biocompatibility tests, need to be performed on the magnetic SIBS material
- ✓ Bonding experiments will be carried out to evaluate the influence of pressure, temperature and time on the bonding strength of the housing element
- ✓ A technique to connect the drainage tube with the housing element will need to be developed

# Acknowledgments

- ✓ Microsystems group at TU/e
- ✓ SEAMS project team



**Maastricht UMC+**



Maastricht University

**Santen**

INNFOCUS

- ✓ This research was conducted under the framework of the Chemelot Institute for Science & Technology (InSciTe)



Chemelot  
**InSciTe**

Thank you for your attention!

